

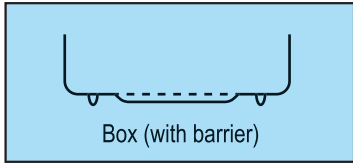


Sterile Container System


3/4 Quarter Size

Bio-Barrier Model

| Specifications Sizes | Lid (with barrier) | Lid (with barrier) | Safety Lid Lid (with barrier) | Safety Lid Lid (with barrier) |
|--|---|---|---|---|
| | Box (Non-Perforated) | Box (with barrier) | Box (Non-Perforated) | Box (with barrier) |
| 465 X 280 X h mm | 11-240 11-241 11-242 11-243 11-244 11-245 h= 115 mm | 11-260 11-261 11-262 11-263 11-264 11-265 h= 135 mm | 11-350 11-351 11-352 11-353 11-354 11-355 h= 128 mm | 11-370 11-371 11-372 11-373 11-374 11-375 h= 150 mm |
| 465 X 280 X h mm | 11-246 11-247 11-248 11-249 11-250 11-251 h= 140 mm | 11-266 11-267 11-268 11-269 11-270 11-271 h= 160 mm | 11-356 11-357 11-358 11-359 11-360 11-361 h= 153 mm | 11-376 11-377 11-378 11-379 11-380 11-381 h= 175 mm |
| 465 X 280 X h mm | 11-252 11-253 11-254 11-255 11-256 11-257 h= 155 mm | 11-272 11-273 11-274 11-275 11-276 11-277 h= 175 mm | 11-362 11-363 11-364 11-365 11-366 11-367 h= 168 mm | 11-382 11-383 11-384 11-385 11-386 11-387 h= 185 mm |
| Ps: All container measurements cover the outer size of the lids. | | | | |
| | | | | |
| | | | | |



| CODE | OUTSIDE | INSIDE | WEIGHT |
|---------------|--------------------|--------------------|---------|
| 11-300 | 475 X 275 X 115 mm | 422 X 260 X 90 mm | 1.95 kg |
| 11-301 | 475 X 275 X 136 mm | 422 X 260 X 111 mm | 2.10 kg |
| 11-302 | 475 X 275 X 157 mm | 422 X 260 X 132 mm | 2.20 kg |

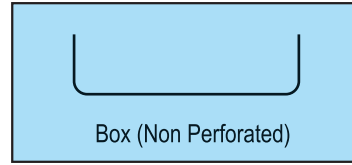

| CODE | OUTSIDE | INSIDE | WEIGHT |
|---------------|--------------------|--------------------|---------|
| 11-305 | 475 X 275 X 95 mm | 422 X 260 X 90 mm | 1.50 kg |
| 11-306 | 475 X 275 X 116 mm | 422 X 260 X 111 mm | 1.65 kg |
| 11-307 | 475 X 275 X 137 mm | 422 X 260 X 132 mm | 1.75 kg |



Lid Barrier System

Weight: 1.285 kg

| | |
|--------|--|
| 11-340 | |
| 11-341 | |
| 11-342 | |
| 11-343 | |
| 11-344 | |
| 11-345 | |

Safety Lids

Weight: 0.94 kg

| | |
|--------|--|
| 11-406 | |
| 11-407 | |
| 11-408 | |
| 11-409 | |
| 11-410 | |
| 11-411 | |






Sterile Container System


1/2 Half Size

Bio-Barrier Model

| Specifications Sizes | Lid (with barrier) Box (Non-Perforated) | Lid (with barrier) Box (with barrier) | Safety Lid Lid (with barrier) Box (Non-Perforated) | Safety Lid Lid (with barrier) Box (with barrier) |
|--|---|---|---|---|
| | 285 X 280 X h mm | 11-500 11-501 11-502 11-503 11-504 11-505 h= 105 mm | 11-530 11-531 11-532 11-533 11-534 11-535 h= 135 mm | 11-640 11-641 11-642 11-643 11-644 11-645 h= 128 mm |
| 285 X 280 X h mm | 11-506 11-507 11-508 11-509 11-510 11-511 h= 140 mm | 11-536 11-537 11-538 11-539 11-540 11-541 h= 160 mm | 11-646 11-647 11-648 11-649 11-650 11-651 h= 153 mm | 11-676 11-677 11-678 11-679 11-680 11-681 h= 170 mm |
| 285 X 280 X h mm | 11-512 11-513 11-514 11-515 11-516 11-517 h= 155 mm | 11-542 11-543 11-544 11-545 11-546 11-547 h= 175 mm | 11-652 11-653 11-654 11-655 11-656 11-657 h= 168 mm | 11-682 11-683 11-684 11-685 11-686 11-687 h= 185 mm |
| 285 X 280 X h mm | 11-518 11-519 11-520 11-521 11-522 11-523 h= 205 mm | 11-548 11-549 11-550 11-551 11-552 11-553 h= 235 mm | 11-658 11-659 11-660 11-661 11-662 11-663 h= 218 mm | 11-688 11-689 11-690 11-691 11-692 11-693 h= 245 mm |
| 285 X 280 X h mm | 11-524 11-525 11-526 11-527 11-528 11-529 h= 265 mm | 11-554 11-555 11-556 11-557 11-558 11-559 h= 295 mm | 11-664 11-665 11-666 11-667 11-668 11-669 h= 278 mm | 11-694 11-695 11-696 11-697 11-698 11-699 h= 305 mm |
| Ps: All container measurements cover the outer size of the lids. | | | | |
| | | | | |

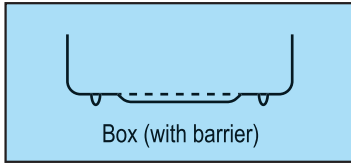

| CODE | OUTSIDE | INSIDE | WEIGHT |
|---------------|--------------------|--------------------|---------|
| 11-600 | 306 X 272 X 115 mm | 253 X 258 X 90 mm | 1.60 kg |
| 11-601 | 306 X 272 X 134 mm | 253 X 258 X 109 mm | 1.70 kg |
| 11-602 | 306 X 272 X 155 mm | 253 X 258 X 130 mm | 1.80 kg |
| 11-603 | 306 X 272 X 204 mm | 253 X 258 X 179 mm | 2.00 kg |
| 11-604 | 306 X 272 X 261 mm | 253 X 258 X 236 mm | 2.45 kg |

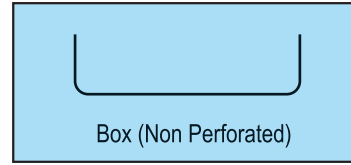

| CODE | OUTSIDE | INSIDE | WEIGHT |
|---------------|--------------------|--------------------|---------|
| 11-605 | 306 X 272 X 95 mm | 253 X 258 X 90 mm | 1.15 kg |
| 11-606 | 306 X 272 X 114 mm | 253 X 258 X 109 mm | 1.25 kg |
| 11-607 | 306 X 272 X 135 mm | 253 X 258 X 130 mm | 1.35 kg |
| 11-608 | 306 X 272 X 184 mm | 253 X 258 X 179 mm | 1.55 kg |
| 11-609 | 306 X 272 X 241 mm | 253 X 258 X 236 mm | 2.00 kg |



Lid Barrier system

Weight: 0.97 kg

| | |
|---------------|--|
| 11-630 | |
| 11-631 | |
| 11-632 | |
| 11-633 | |
| 11-634 | |
| 11-635 | |

Safety Lids

Weight: 0.66 kg

| | |
|---------------|--|
| 11-412 | |
| 11-413 | |
| 11-414 | |
| 11-415 | |
| 11-416 | |
| 11-417 | |



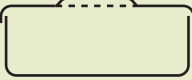

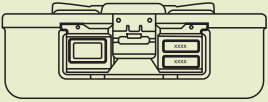












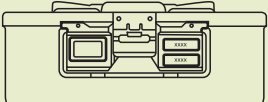















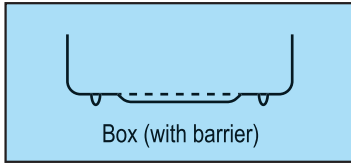
Sterile Container System

.


XL Size

Bio-Barrier Model

| <p>Specifications</p> <p>Sizes</p> | <p>Lid (with barrier)</p>  <p>Box (Non-Perforated)</p> | <p>Lid (with barrier)</p>  <p>Box (with barrier)</p> |
|---|--|--|
| <p>590 X 370 X h mm</p>  | <p>11-711 </p> <p>11-712 </p> <p>11-713 </p> <p>11-714 </p> <p>11-715 </p> <p>11-716 </p> <p style="text-align: right;">h = 125 mm</p> | <p>11-741 </p> <p>11-742 </p> <p>11-743 </p> <p>11-744 </p> <p>11-745 </p> <p>11-746 </p> <p style="text-align: right;">h = 135 mm</p> |
| <p>590 X 370 X h mm</p>  | <p>11-717 </p> <p>11-718 </p> <p>11-719 </p> <p>11-720 </p> <p>11-721 </p> <p>11-722 </p> <p style="text-align: right;">h = 185 mm</p> | <p>11-747 </p> <p>11-748 </p> <p>11-749 </p> <p>11-750 </p> <p>11-751 </p> <p>11-752 </p> <p style="text-align: right;">h = 195 mm</p> |
| <p>Ps: All container measurements cover the outer size of the lids.</p> |  | |



| CODE | OUTSIDE | INSIDE | WEIGHT |
|---------------|--------------------|--------------------|----------|
| 11-700 | 595 X 365 X 125 mm | 530 X 350 X 100 mm | 2.620 kg |
| 11-701 | 595 X 365 X 187 mm | 530 X 350 X 162 mm | 3.440 kg |

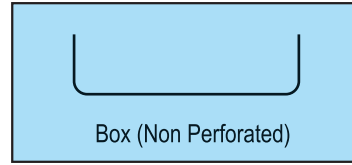

| CODE | OUTSIDE | INSIDE | WEIGHT |
|---------------|--------------------|--------------------|---------|
| 11-702 | 595 X 365 X 105 mm | 530 X 350 X 100 mm | 2.20 kg |
| 11-703 | 595 X 365 X 167 mm | 530 X 350 X 167 mm | 3.00 kg |




Lid Barrier system


Weight: 2.037 kg

| | |
|---------------|--|
| 11-771 | |
| 11-772 | |
| 11-773 | |
| 11-774 | |
| 11-775 | |
| 11-776 | |




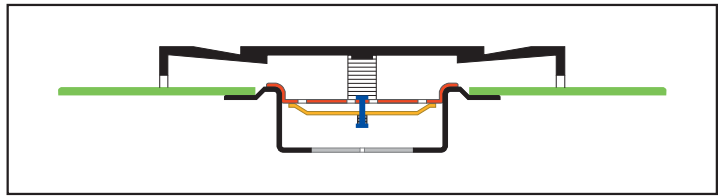

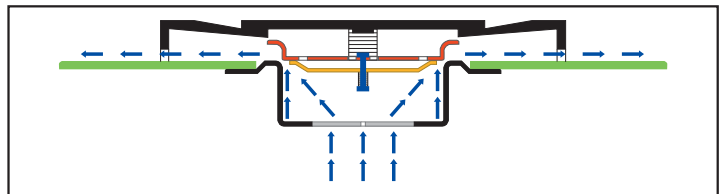
1. Before sterilization the valves are in relax position.



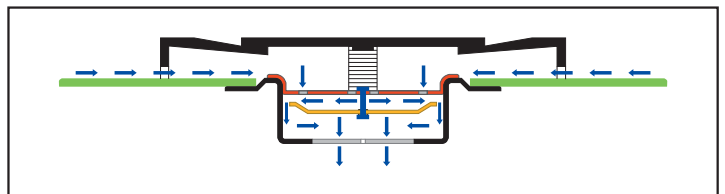
2. During the sterilization by effect of vacuum the valves are opened completely and the dry air goes out of the container.



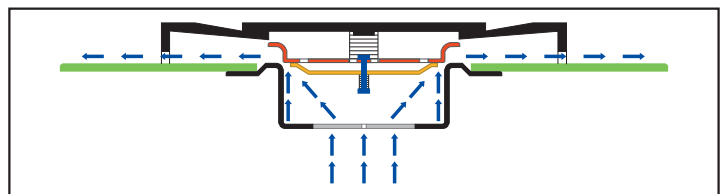3. At the end the sterilization cycle the valves are opened in the reverse way by effect of vacuum and the sterilization cycle is completed.



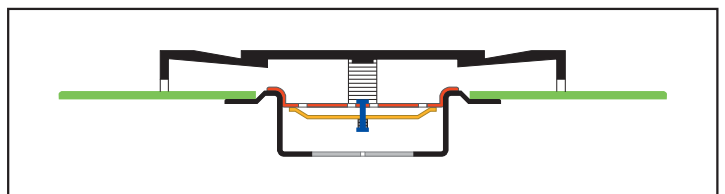4. At the end of the sterilization cycle the sterilizer equates the pressure level and the valves of the Barrier Model Container are closed.



5. After sterilization the valves are in relax position.



Sterile Container System

Dental Mini & Special Model

| Specifications Sizes | Lid (Perforated) | Lid (Perforated) | Lid (Non Perforated) | | |
|-----------------------------|--|--|--|------------------|----------------------|
| | | | | | |
| | Box (Non-Perforated) | Box (Perforated) | Box (Non-Perforated) | Box (Perforated) | Box (Non-Perforated) |
| 320 X 190 X 50 mm | 10-1800 10-1801 10-1802 10-1803 10-1804 10-1805 | 10-1820 10-1821 10-1822 10-1823 10-1824 10-1825 | 10-1840 10-1841 10-1842 10-1843 10-1844 10-1845 | 10-1860 | 10-1865 |
| 320 X 190 X 65 mm | 10-1806 10-1807 10-1808 10-1809 10-1810 10-1811 | 10-1826 10-1827 10-1828 10-1829 10-1830 10-1831 | 10-1846 10-1847 10-1848 10-1849 10-1850 10-1851 | 10-1861 | 10-1866 |
| 320 X 190 X 130 mm | 10-1812 10-1813 10-1814 10-1815 10-1816 10-1817 | 10-1832 10-1833 10-1834 10-1835 10-1836 10-1837 | 10-1852 10-1853 10-1854 10-1855 10-1856 10-1857 | 10-1862 | 10-1867 |

Ps: All container measurements cover the outer size of the lids.



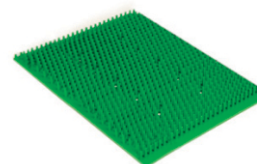
13-1070 274 x 172 x 43 mm
WIRE BASKET



13-611 274 x 172 x 43 mm
WIRE BASKET WITH LID.



13-609 273 x 180 x 13 mm
STORAGE FRAME for dental model




12-529 260 x 180 mm
SILICONE MAT



13-610 277 x 175 mm
DRAPE RETAINER


Perforated Bottoms



| |
|------------------------------|
| outside 320 x 186 x 38 mm |
| inside 276 x 185 x 23 mm |
| weight 0,45 kg |


10-1860

Solid Bottoms




| |
|------------------------------|
| outside 320 x 186 x 38 mm |
| inside 276 x 185 x 36 mm |
| weight 0,35 kg |

10-1865




| |
|------------------------------|
| outside 320 x 186 x 55 mm |
| inside 276 x 185 x 40 mm |
| weight 0,50 kg |

10-1861




| |
|------------------------------|
| outside 320 x 186 x 55 mm |
| inside 276 x 185 x 43 mm |
| weight 0,40 kg |

10-1866



| |
|-------------------------------|
| outside 320 x 186 x 123 mm |
| inside 276 x 185 x 108 mm |
| weight 0,70 kg |


10-1862



| |
|-------------------------------|
| outside 320 x 186 x 123 mm |
| inside 276 x 185 x 121 mm |
| weight 0,65 kg |

10-1867

Perforated Lids



| | |
|---------|--|
| 10-1880 | |
| 10-1881 | |
| 10-1882 | |
| 10-1883 | |
| 10-1884 | |
| 10-1885 | |

Weight: 0.50 kg




13-615
Single use paper filter 100 pcs.
With indicator No marking.



12-743
Teflon filter
Number of sterilization cycle > 1200


Bio-Barrier Lids



| | |
|---------|--|
| 11-1880 | |
| 11-1881 | |
| 11-1882 | |
| 11-1883 | |
| 11-1884 | |
| 11-1885 | |


Weight: 0.90 kg




Removal Lids




| | |
|---------|--|
| 10-1870 | |
| 10-1871 | |
| 10-1872 | |
| 10-1873 | |
| 10-1874 | |
| 10-1875 | |

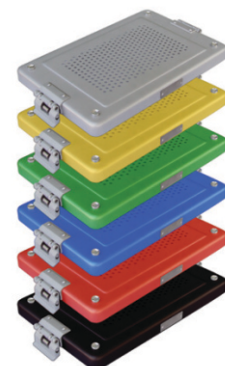Weight: 0.60 kg




| Specifications Sizes | Lid (Perforated) | Lid (Perforated) | Lid (Non Perforated) | | |
|---|--|--|--|----------------------------------|----------------------|
| | Box (Non-Perforated) | Box (Perforated) | Box (Non-Perforated) | Box (Perforated) | Box (Non-Perforated) |
| 315 X 135 X 40 mm  | 10-1900 10-1901 10-1902 10-1903 10-1904 10-1905 | 10-1920 10-1921 10-1922 10-1923 10-1924 10-1925 | 10-1940 10-1941 10-1942 10-1943 10-1944 10-1945 | 10-1960 10-1965 | |
| 315 X 135 X 70 mm  | 10-1906 10-1907 10-1908 10-1909 10-1910 10-1911 | 10-1926 10-1927 10-1928 10-1929 10-1930 10-1931 | 10-1946 10-1947 10-1948 10-1949 10-1950 10-1951 | 10-1961 10-1966 | |
| 315 X 135 X 100 mm  | 10-1912 10-1913 10-1914 10-1915 10-1916 10-1917 | 10-1932 10-1933 10-1934 10-1935 10-1936 10-1937 | 10-1952 10-1953 10-1954 10-1955 10-1956 10-1957 | 10-1962 10-1967 | |

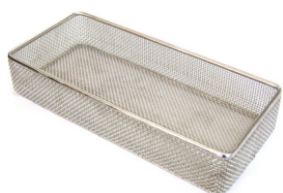Ps: All container measurements cover the outer size of the lids.



13-1073 267 x 125 x 50 mm
WIRE BASKET

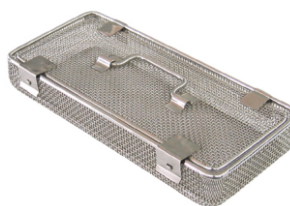

12-993 267 x 125 x 50 mm
WIRE BASKET WITH LID



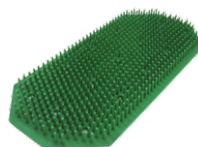
12-990 270 x 125 mm
DRAPE RETAINER



12-991
Portable carrier for small containers. Wire suspenders adjustable as required. Adjusting key of suspender is included. Width: 325 mm, Height: 328 mm, Depth: 145 mm




12-992 270 x 125 mm
WIRE BASE



12-520 275 x 125 mm
SILICONE MAT


Perforated Bottoms



| |
|------------------------------|
| outside 315 x 132 x 32 mm |
| inside 273 x 128 x 24 mm |
| weight 0,30 kg |


10-1960

Solid Bottoms




| |
|------------------------------|
| outside 296 x 132 x 32 mm |
| inside 273 x 128 x 30 mm |
| weight 0,20 kg |

10-1965




| |
|------------------------------|
| outside 315 x 132 x 58 mm |
| inside 273 x 128 x 50 mm |
| weight 0,40 kg |

10-1961




| |
|------------------------------|
| outside 296 x 132 x 58 mm |
| inside 273 x 128 x 56 mm |
| weight 0,35 kg |

10-1966



| |
|------------------------------|
| outside 315 x 132 x 92 mm |
| inside 273 x 128 x 84 mm |
| weight 0,50 kg |

10-1962



| |
|------------------------------|
| outside 296 x 132 x 92 mm |
| inside 273 x 128 x 90 mm |
| weight 0,45 kg |

10-1967




12-761
Single use paper filter 100 pcs.
With indicator. No marking.



12-744
Teflon filter
Number of sterilization scycle > 1200


Perforated Lid



| | |
|---------|--|
| 10-1970 | |
| 10-1971 | |
| 10-1972 | |
| 10-1973 | |
| 10-1974 | |
| 10-1975 | |

Weight: 0.35 kg

Removal Lids



| | |
|---------|--|
| 10-1980 | |
| 10-1981 | |
| 10-1982 | |
| 10-1983 | |
| 10-1984 | |
| 10-1985 | |

Weight: 0.30 kg



IMPLANT Model Containers

| Specifications / Sizes | Lid (Perforated) | Lid (Perforated) | Lid (Non Perforated) | NON FILTER |
|---------------------------|--|--|--|------------|
| | Box (Non-Perforated) | Box (Perforated) | Box (Non-Perforated) | |
| 300 X 149 X 65 mm | 09-010 09-011 09-012 09-013 09-014 09-015 | 09-022 09-023 09-024 09-025 09-026 09-027 | 09-034 09-035 09-036 09-037 09-038 09-039 | |
| 500 X 169 X 75 mm | 09-016 09-017 09-018 09-019 09-020 09-021 | 09-028 09-029 09-030 09-031 09-032 09-033 | 09-040 09-041 09-042 09-043 09-044 09-045 | |


Ps: All container measurements cover the outer size of the lids.

MINI IMPLANT Model Containers

| Specifications / Sizes | Lid (Perforated) | Lid (Perforated) | Lid (Non Perforated) | NON FILTER |
|---------------------------|--|--|--|------------|
| | Box (Non-Perforated) | Box (Perforated) | Box (Non-Perforated) | |
| 160 X 70 X 60 mm | 09-200 09-201 09-202 09-203 09-204 09-205 | 09-212 09-213 09-214 09-215 09-216 09-217 | 09-224 09-225 09-226 09-227 09-228 09-229 | |
| 200 X 145 X 60 mm | 09-206 09-207 09-208 09-209 09-210 09-211 | 09-218 09-219 09-220 09-221 09-222 09-223 | 09-230 09-231 09-232 09-232 09-234 09-235 | |

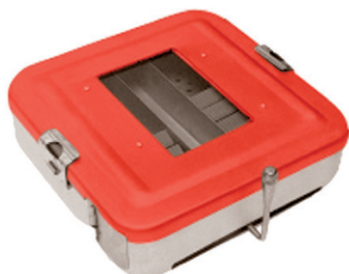Ps: All container measurements cover the outer size of the lids.



**BRUSH Model
Containers**

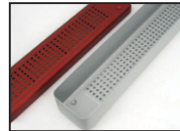
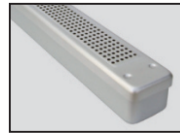


- 12-044 Black
- 12-045 Red
- 12-046 Blue
- 12-047 Green
- 12-048 Yellow
- 12-049 Grey

Container for brush 280 x 285 x 85 mm with wall mounting fixtures, Maximum 20 each brushes can be stored container to take the brush, please move the handle to right or left direction under the container.



non filter


| Specifications / Sizes | Lid (Perforated) | Lid (Perforated) |
|--|--|--|
| | Box (Non-Perforated) | Box (Perforated) |
| 250 X 60 X 50 mm  | 09-060 09-061 09-062 09-063 09-064 09-065 | 09-072 09-073 09-074 09-075 09-076 09-077 |
| 450 X 70 X 70 mm  | 09-066 09-067 09-068 09-069 09-070 09-071 | 09-078 09-079 09-080 09-081 09-082 09-083 |



SCOPY Model Containers



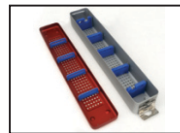Ps: All container measurements cover the outer size of the lids.

| Specifications / Sizes | Lid (Perforated) | Lid (Perforated) |
|---|--|--|
| | Box (Non-Perforated) | Box (Perforated) |
| 250 X 60 X 50 mm  | 09-100 09-101 09-102 09-103 09-104 09-105 | 09-115 09-116 09-117 09-118 09-119 09-120 |
| 450 X 70 X 70 mm  | 09-106 09-107 09-108 09-109 09-110 09-111 | 09-121 09-122 09-123 09-124 09-125 09-126 |



SCOPY Model Containers
(with silicon stand, lock, label)

Ps: All container measurements cover the outer size of the lids.

| Specifications Sizes | Lid (Perforated) | Lid (Perforated) | Lid (Non Perforated) | | |
|--|--|--|--|------------------|----------------------|
| | Box (Non-Perforated) | Box (Perforated) | Box (Non-Perforated) | Box (Perforated) | Box (Non-Perforated) |
| 285 X 280 X 55 mm  | 10-1710 10-1711 10-1712 10-1713 10-1714 10-1715 | 10-1730 10-1731 10-1732 10-1733 10-1734 10-1735 | 10-1750 10-1751 10-1752 10-1753 10-1754 10-1755 | 10-1770 | 10-1775 |
| 285 X 280 X 85 mm  | 10-1716 10-1717 10-1718 10-1719 10-1720 10-1721 | 10-1736 10-1737 10-1738 10-1739 10-1740 10-1741 | 10-1756 10-1757 10-1758 10-1759 10-1760 10-1761 | 10-1771 | 10-1776 |

Ps: All container measurements cover the outer size of the lids.



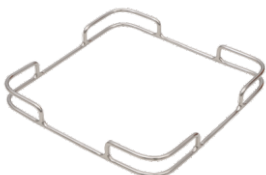
12-490 255 x 250 x 30 mm
12-491 255 x 250 x 50 mm
 WIRE BASKET




12-482 255 x 255 x 30 mm
 WIRE BASE



12-480 255 x 255 mm
 DRAPE RETAINER



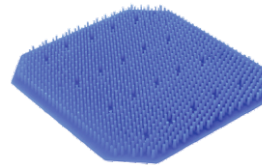12-481 255 x 255 x 30 mm
 DRAPE RETAINER



12-521 255 x 255 mm
 SILICON MAT




12-740
 Single use paper filter 100 pcs.



12-742
 Teflon Filter
 Number of sterilization cycle > 1200


Perforated Bottoms



| |
|------------------------------|
| outside 209 x 275 x 40 mm |
| inside 264 x 270 x 19 mm |
| weight 0,70 kg |


10-1770

Solid Bottoms



| |
|-----------------------------|
| outside 290x 275 x 40 mm |
| inside 264 x 270 x 35 mm |
| weight 0,60 kg |

10-1775



| |
|------------------------------|
| outside 290 x 275 x 72 mm |
| inside 264 x 270 x 52 mm |
| weight 0,90 kg |


10-1771



| |
|------------------------------|
| outside 290 x 275 x 72 mm |
| inside 264 x 270 x 66 mm |
| weight 0,70 kg |

10-1776


Perforated Lids



| | |
|---------|--|
| 10-1780 | |
| 10-1781 | |
| 10-1782 | |
| 10-1783 | |
| 10-1784 | |
| 10-1785 | |

Weight: 0.65 kg


Bio-Barrier Lids



| | |
|---------|--|
| 11-1780 | |
| 11-1781 | |
| 11-1782 | |
| 11-1783 | |
| 11-1784 | |
| 11-1785 | |

Weight: 0.85 kg


Safety Lids



| | |
|---------|--|
| 10-1418 | |
| 10-1419 | |
| 10-1420 | |
| 10-1421 | |
| 10-1422 | |
| 10-1423 | |

Weight: 0.55 kg


Bio-Barrier Safety Lids



| | |
|---------|--|
| 11-1418 | |
| 11-1419 | |
| 11-1420 | |
| 11-1421 | |
| 11-1422 | |
| 11-1423 | |


Weight: 0.60 kg

Removal Lids



| | |
|---------|--|
| 10-1790 | |
| 10-1791 | |
| 10-1792 | |
| 10-1793 | |
| 10-1794 | |
| 10-1795 | |

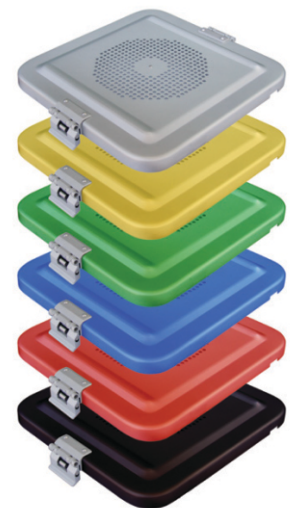Weight: 0.50 kg



Container System & Baskets

Classic & Bio-Barrier Model



Surgical innovation Sytems

From Access to Repair

Address:

Veko Giz Plaza, Meydan Sokak No.3/43, Kat.13 Maslak, Istanbul,Turkey

Phone: +90 212 705 36 16, Fax: +90 212 705 36 36

info@thempson.com

www.thempson.com